



**silfradent**<sup>®</sup>

## SURGYBONE TIPS

A new generation in ultrasound surgery



**SILFRADENT**, with its 30 years of experience, manufactures products in compliance with the European regulations.

# SURGYBONE TIPS

## SB P0100 ✓ CUTTING TIPS

Scalpel angled at 90°. Ideally suited for sinus lift and bone harvesting.



✓ BONE SCRAPER



- PWR ULT 30
- VIBRA 100

## SB P0200 ✓ CUTTING TIPS

Non-traumatic diamond tip. Ideally suited when opening the window.



✓ SINUS LIFT WINDOW



- PWR ULT 30
- VIBRA 100

## SB P0301 ✓ CUTTING TIPS

Round shaped cutting scalpel, angled at 120°. Ideally suited for periodontics and bone harvesting.



✓ BONE RECOVERY



- PWR ULT 30
- VIBRA 100

✓ PERIODONTICS




- PWR ULT 35
- VIBRA 100

## SB P0310 ✓ NON-CUTTING TIPS

To lift the sinus membrane, angled at 90°, flat circular.



✓ SINUS LIFT MEMBRANE

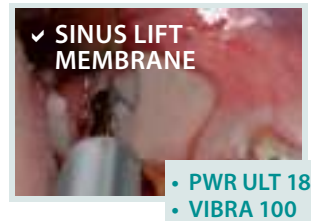

- PWR ULT 18
- VIBRA 100

## SB P0321 ✓ NON-CUTTING TIPS

To lift the sinus membrane, angled at 120°, flat circular.



✓ SINUS LIFT MEMBRANE



- PWR ULT 20
- VIBRA 100

## SB P0400 ✓ NON-CUTTING TIPS

To separate the bone from the maxillary sinus.




✓ SINUS LIFT MEMBRANE



- PWR ULT 20
- VIBRA 100

### SB P0500 ✓ CUTTING TIPS

Cutting tip with scalpel top and ridged sides.  
Ideally suited when harvesting bone.

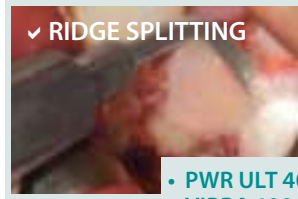

### SB P0600 ✓ CUTTING TIPS

Concaved cutting tip.  
Ideally suited when harvesting bone from the mandible.



### SB P0610 ✓ CUTTING TIPS

120° sharp lance tip for extraction and removal of teeth.




### SB P0700 ✓ CUTTING TIPS

Fast and effective cutting tip.



### SB P0710 ✓ CUTTING TIPS

Fast and effective cutting tip.  
Angled towards the right.

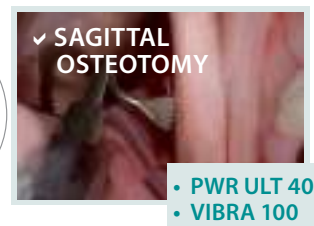

### SB P0720 ✓ CUTTING TIPS

Fast and effective cutting tip.  
Angled towards the left.




**SB P0750** ✓ CUTTING TIPS

Fast and effective cutting tip.




**SB P0800** ✓ SMOOTHING TIPS

Non-traumatic ballend diamond tip. Ideally suited when working close to nerve endings and blood vessels.

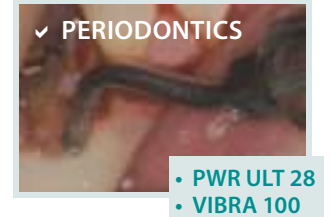



**SB P0805** ✓ SMOOTHING TIPS

Non-traumatic diamond tip. Ideally suited on thin bone and cartilage.

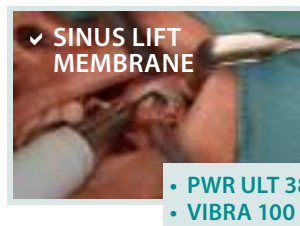


**SB P0830** ✓ CUTTING TIPS

Non-traumatic microsharpening tip. Ideally suited to prepare and compact the implant site and in case of extraction implant.

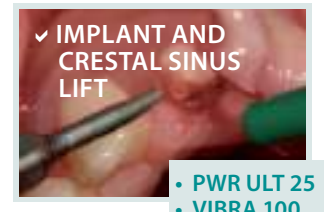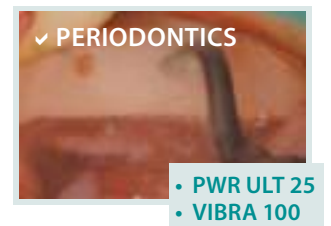



**SB P0900** ✓ CUTTING TIPS

Non-traumatic polished cone shaped tip. Ideally suited for periodontal surgery.

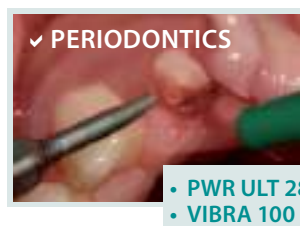

**SB P0940** ✓ CUTTING TIPS

Non-traumatic microsharpening tip. Ideally suited to prepare and compact the implant site and in case of post extraction implant.

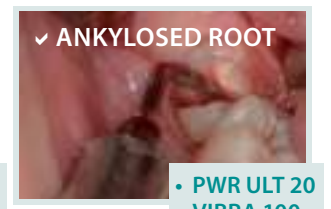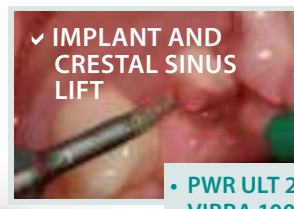

# SURGYBONE TIPS

## KIT PEC // COMPLETE WITH 6 TIPS



## KIT USE HYDRO // COMPLETE WITH 6 TIPS



## KIT MAXILLO



## KIT SPYNE SURGERY



- ✓ PWR ULT: ULTRASOUND POWER
- ✓ VIBRA: CUTTING POWER

This brochure contains the average standard parameters provided by Silfradent which can be adjusted according to the several variables to be considered in surgery. We invite you to join us in our training courses. Our tutors will help you understand why the Surgybone is presently the highest performance ultrasound machine.

## SURGYBONE

### A METHODOLOGY FOR OSTEOIMPLANT SURGERY

Technology based on ultrasound frequencies that permits highly precise osteotomies, with a constant control of both incision depth and length.



#### AREA OF USE

- ✓ Sinus Lift with crestal approach
- ✓ Extraction
- ✓ Apicectomy
- ✓ Cystectomy
- ✓ Bone Harvesting
- ✓ Ridge Expansion
- ✓ Sinus Lift
- ✓ Osteoplasmy, Osteotomy
- ✓ Implant site preparation
- ✓ Tooth preparation and root planning
- ✓ Orthodontic Surgery

**OUR RESEARCH  
IS JUST THE RESULT  
OF YOUR  
SUGGESTIONS**

A product can be developed and improved just through the experience derived from practice. Silfradent Medical Line believes in the satisfaction of user's needs. For this reason and in order to guarantee a better service, we listen to your suggestions, that you can send directly to [info@silfradent.com](mailto:info@silfradent.com).



Equipment for **Medical Line**



# ODONTES

DENTAL AND MEDICAL SOLUTIONS

Distribuito da:  
ODONTES SRL - Via Don Sturzo,6 - 20016 Pero (MI)  
T. 02.38.10.30.89 - [info@odontes.it](mailto:info@odontes.it)