

STRIP CROWNS




ODONTES  
DENTAL AND MEDICAL SOLUTIONS

Distribuito da:  
ODONTES SRL - Via Don Sturzo,6 - 20016 Pero (MI)  
T. 02.38.10.30.89 - info@odontes.it

## frasaco-strip crowns

- \* renowned matrices for fast and clean temporary restorations
- \* elastic, thin and easy to cut, packed antistatically
- \* pedodontal strip crown available for restoration of badly damaged deciduous teeth
- \* wide range covering all required sizes
- \* assortments for proper storage and choice, refill-packages available

Save time and money with frasaco strip crowns, achieving optimum efficiency in no time. Even for multiple layer technique, strip crowns can be used as matrix to shape the final layer.

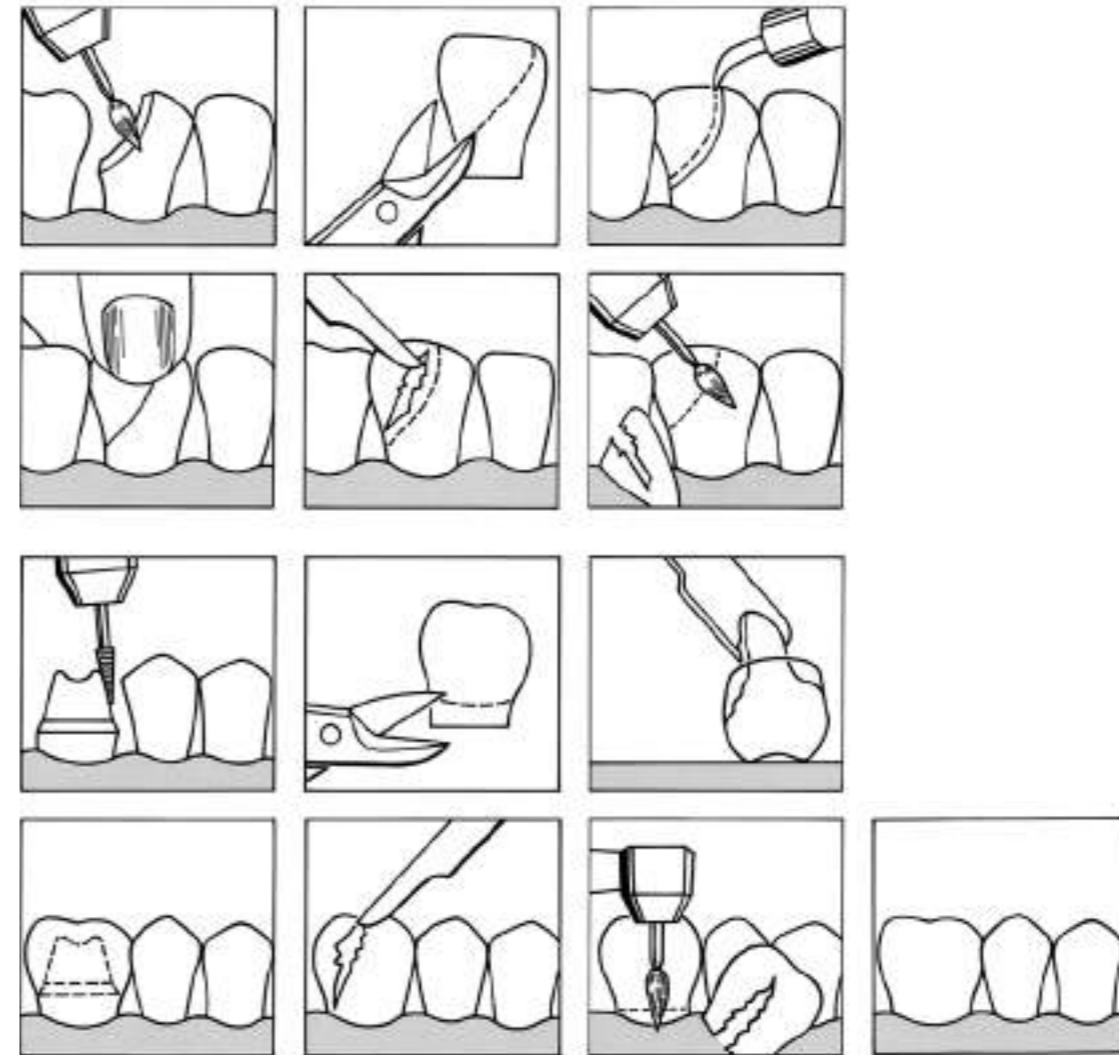

frasaco strip crowns are elastic, thin and easy to cut celluloid crowns. In the transparent execution they are light-transmissive for polymerization light. They are used as moulds for creating temporary restorations on preparations and as matrices for restoration of corners and edges. They are placed without problem using any conventional material for temporary restorations. A strip crown in the required size and shape is adjusted to the cervical area and the gingival margin with a crown scissor. Please request manufacturers' addresses for products we have tested for this application.

Prior to placement, ensure that strip crown and preparation are clean and dry.

Fill the strip crown with temporary filling material and place it. Strip crowns can be used as matrices when restoring corners and edges – even when using composites. When doing so, cut out that section of the strip crown corresponding to the part of the tooth being restored, allowing for the strip crown to overlap the margins of the cavity. Fill the matrix with the definite filling material and press it into place. Remove the matrix after the filling material has set. Should for some reason strip crowns remain for a few days on a tooth, please make sure that their gingival edge is no deeper than 1mm inside the sulcus to avoid possible traumatization of the marginal periodontium. Under no circumstances should strip crowns remain for longer than 30 days in the mouth.

# strip crowns – forms and sizes

central incisors  
– right upper



central incisors  
left upper –

lateral incisors  
– right upper



lateral incisors  
left upper –

canines  
– right upper



canines  
left upper –

first premolars  
– right upper



first premolars  
left upper –

first molars  
– right upper



first molars  
left upper –

first molar  
– right lower



first molars  
left lower –

first premolars  
– right lower



first premolars  
left lower –

canines  
– right lower



canines  
left lower –

lateral incisors  
– right lower



lateral incisors  
left lower –

central incisors  
– right lower



central incisors  
left lower –

## transparent strip crowns



### ST-5K

5 transparent strip crowns of any one size  
refill package

please request our order form for strip crowns



### ST-FK174

anterior assortment  
174 transparent strip crowns  
plastic box with transparent cover  
3 each of every size in 58 compartments


Upper 110-115	Lower 311-313
Upper 120-125	Lower 320-322
Upper 130-135	Lower 330-334
Upper 210-215	Lower 411-413
Upper 220-225	Lower 420-422
Upper 230-235	Lower 430-434



### ST-BK96

posterior assortment  
96 transparent strip crowns  
plastic box with transparent cover  
3 each of every size in 32 compartments

Upper 141-144	Lower 341-344
Upper 160-163	Lower 360-363
Upper 241-244	Lower 441-444
Upper 260-263	Lower 460-463

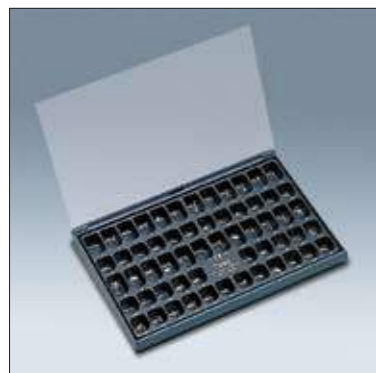

### ST-FBK 58

compact assortment (current anteriors and posteriors)  
58 transparent strip crowns  
plastic box with transparent cover  
1 strip crown of each one size in 58 compartments

Upper 111-116	Lower 311-312
Upper 120-124	Lower 321-322
Upper 131, 133	Lower 330-331
Upper 141-143	Lower 341-342
Upper 161-162	Lower 360-362
Upper 211-216	Lower 411-412
Upper 220-224	Lower 421-422
Upper 231, 233	Lower 430-431
Upper 241-243	Lower 441-442
Upper 261-262	Lower 460-462

## coloured strip crowns

### ST-5E

5 coloured strip crowns of any one size  
refill package


please request our order form for strip crowns



### ST-FE174

anterior assortment  
174 coloured strip crowns  
plastic box with transparent cover  
3 each of every size in 58 compartments


Upper 110-115	Lower 311-313
Upper 120-125	Lower 320-322
Upper 130-135	Lower 330-334
Upper 210-215	Lower 411-413
Upper 220-225	Lower 420-422
Upper 230-235	Lower 430-434

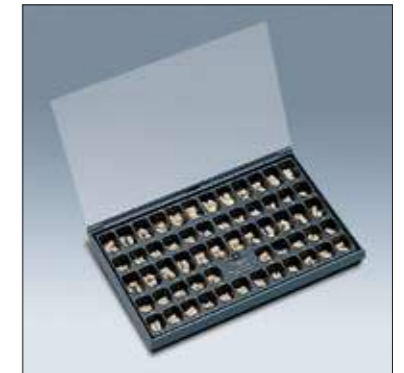

### ST-BE96

posterior assortment  
96 coloured strip crowns  
plastic box with transparent cover  
3 each of every size in 32 compartments

Upper 141-144	Lower 341-344
Upper 160-163	Lower 360-363
Upper 241-244	Lower 441-444
Upper 260-263	Lower 460-463

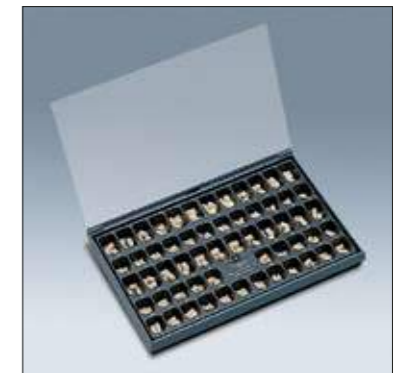

### ST-FBE 58

compact assortment (current anteriors and posteriors)  
58 coloured strip crowns  
plastic box with transparent cover  
1 strip crown of each one size in 58 compartments

Upper 111-116	Lower 311-312
Upper 120-124	Lower 321-322
Upper 131, 133	Lower 330-331
Upper 141-143	Lower 341-342
Upper 161-162	Lower 360-362
Upper 211-216	Lower 411-412
Upper 220-224	Lower 421-422
Upper 231, 233	Lower 430-431
Upper 241-243	Lower 441-442
Upper 261-262	Lower 460-462



## pedodontal strip crowns

### Pedodontal strip crowns

Placing prefabricated frasaco pedodontal strip crowns is an excellent technique for retaining badly damaged deciduous teeth.

In earlier times, deciduous teeth with carious defects were extracted once these exceeded a certain stage. Nowadays, they are left in the mouth as a spacer as it is now known that the deciduous teeth are important for the unimpeded development of children's jaws. The advances made in conservative dentistry have changed the manner in which this field is approached. Reference should also be made to aesthetic and phonetic aspects.

We have kept pace with these changes by supplying our pedodontal strip crowns in several sizes. Allergies, which can be caused by other crown materials, are virtually eliminated.

Please request manufacturers' addresses for crown scissors which we have tested for this application.

#### central incisors



#### lateral incisors



#### canines



frasaco



**ST-K5K** 5 transparent strip crowns of any one size, refill package

**ST-K5E** 5 coloured strip crowns of any one size, refill package

please request our order form for strip crowns



**ST-KK72** pedodontal strip crown assortment  
72 transparent strip crowns for upper anteriors  
plastic box with transparent cover  
3 strip crowns of each one size in 24 compartments  
Upper 511-514, Upper 521-524  
Upper 531-534, Upper 611-614  
Upper 621-624, Upper 631-634

**ST-KE72** pedodontal strip crown assortment  
72 coloured strip crowns for upper anteriors, plastic box with transparent cover  
3 strip crowns of each one size in 24 compartments